
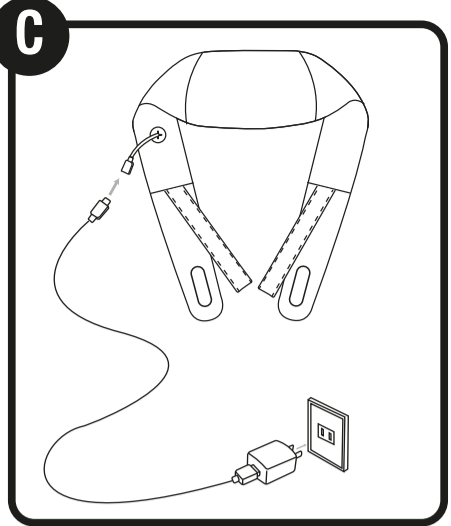
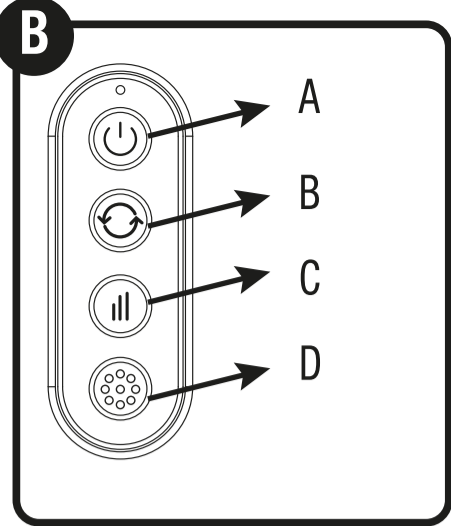
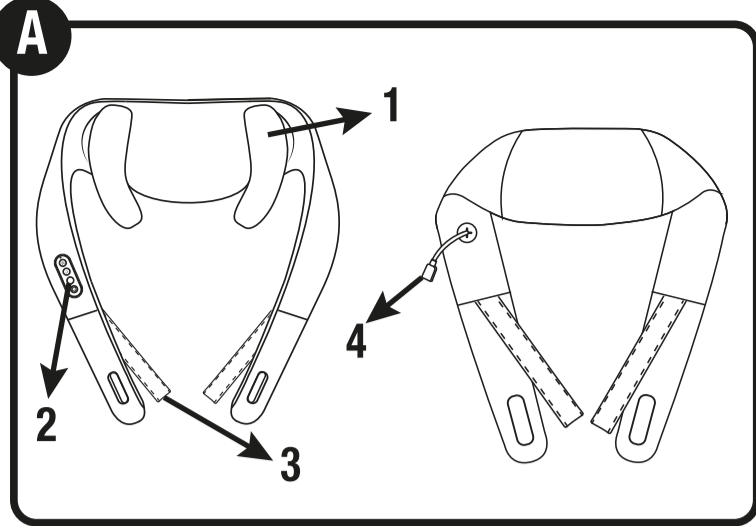
WELLPRIME®

Total Neck Massager

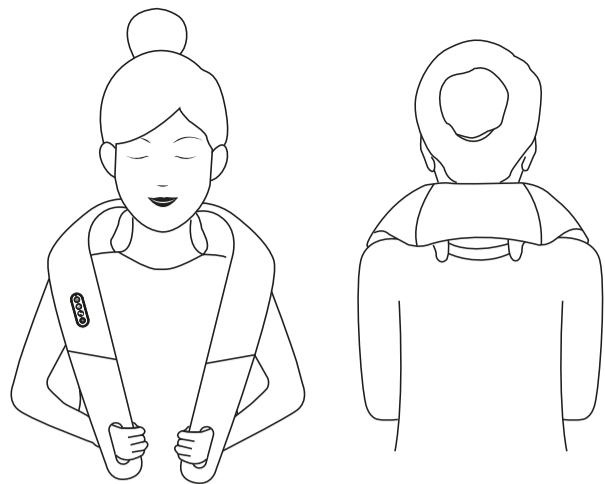ENGLISH

Thank you for purchasing Wellprime® Total Neck Massager, the ergonomically designed shoulder and neck massager that adapts to your neck and keeps your head in a proper position. Its massage heads simulate the feeling of a real hand massage.

Please read the instructions for use and all safety instructions carefully before use. Keep the manual in an accessible place for future use.

SAFETY WARNINGS

1. Do not place the Wellprime® Total Neck Massager in the bathroom or in a damp place.
2. Always keep the Wellprime® Total Neck Massager dry. Do not store it in damp places or near containers of water into which it may fall.
3. It is not waterproof; if the unit should fall into water, do not touch or reuse it. Do not submerge the unit under water.
4. Children under 3 years of age should be kept away unless continuously supervised.
5. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children must not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.

6. Do not leave the massager exposed to direct sunlight or high temperatures for a long period of time.

7. This product has a built-in rechargeable battery. Please dispose of this product in accordance with current environmental protection guidelines.

8- Never charge the battery if the temperature is below 0°C. Use, charge and store the device at a temperature between 5°C and 35°C.

9- Avoid charging the battery while exposed to direct sunlight or near a heat source.

10. Never charge the battery uninterruptedly for more than 48 hours.

11- To recharge the battery, use only the detachable power supply unit supplied with this appliance.

12- The appliance may only be used with the supplied power supply unit.

13- Keep the removed rechargeable battery out of the reach of children to prevent accidental ingestion. In case of accidental ingestion, seek medical advice immediately.

14- In case of battery fluid leakage, do not touch the fluid with unprotected hands and handle it as follows:

If battery liquid enters the eyes, there is a danger of blindness, do not rub the eyes directly with your hands, rinse them with clean water immediately and seek medical attention. If battery fluid comes into contact with the body or clothing, it may cause inflammation or injury to the skin -- rinse with clean water. If you continue to feel unwell, seek medical attention immediately.

15- This product is not a toy.
16. Do not use this product if it is not working properly and the massage device is damaged or broken.

17. In case of severe or intense pain during the massage, stop the unit immediately and seek medical advice.

18. Do not use this product directly on areas that are swollen, inflamed or with skin rashes.
19. Do not charge the product when the power supply is unstable.

20. Do not touch the power supply or plug with wet hands to avoid malfunction or electric shock.

21- Do not use the product lying down or half lying down.

22- Do not use the Wellprime® Total Neck Massager if you have heart problems or heart disease.

23- Patients with severe cervical problems and/or skeletal hyperostosis must not use the product.

24- Patients with spinal problems should not use the product.

25- People with severe blood circulation problems in the shoulder and neck should not use the product.

26- Do not disassemble, repair or modify this product without authorization. Failure to do so may result in malfunction or electric shock.

- 2- COMPONENTS
FIGURE A
1- Massage head
2- Control buttons/ Indicator
3- Strap
4- Type-C Charging Port

- FIGURE B
Control Buttons / Indicator
A- On/Off
B- Mode
C- Speed
D- Hot compression

3-MODE OF USE

1. Power on: Press and hold the button for 2 seconds to turn on the Wellprime® Total Neck Massager.
2. Rotation mode selection: the massage heads rotate in two directions. Briefly press the button to change the direction of the massage.
3. Heat function: By pressing the heat function button the Wellprime® Total Neck Massager will emit heat. Press once for high temperature, press the button briefly again to switch to a lower temperature, and then press briefly one more time to deactivate the heat function -- cycle in this way to repeat the process.
4. Speed function: Wellprime® Total Neck Massager has two three; briefly press the button to switch high, medium or low speed.
5. Switch off
Automatic shut-off: The product stops working after a total of 10 minutes of use.
Manual shut-off: When you stop using the product, press and hold the power button for 2 seconds to turn it off.

NOTE: Wellprime® Total Neck Massager has a strap to adjust the position of the Wellprime® Total Neck Massager.

4- BATTERY CHARGING

The massager is powered by a built-in lithium battery. If the battery is about to or has run out of power, charge the massager before use.

Do not use the product while charging and be sure to plug in the charging cable first.
Insert the plug of the charging cable into the charging interface. Charging status: the red light flashes during charging, and is always on when fully charged.

Wellprime® Total Neck Massager takes about 2 hours to charge. This product includes accessories that are only for use with this machine. Other accessories on the market are not to be used as substitute for the accessories that come with this product.
NOTE: The USB cable is included in the set. The mains plug-in adapter is not included. (See figure C)

5- Cleaning and Maintenance Cover:
After each use: turn off the massager, and store it in a dry and safe place.
Wellprime® Total Neck Massager cannot be washed or immersed in water.

Clean the outer surface of the massager with a damp cloth. Do not rinse the product directly with water.
CAUTION: Never use detergents to clean the massager, as they may damage the operating parts.

TECHNICAL SPECIFICATIONS

Power : 26 w
Input: 5V=1A
Battery: 2000 mah
Charging voltage: 5V

Composition:
Fabric : PU; 90% Polyester+10% spandex
Padding: Viscoelastic foam

QUALITY WARRANTY:

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.
This warranty does not cover damage resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

Australia & New Zealand only
This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

Disposal of old electric appliances
The European directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream.

Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.
Consumers should contact their local authority or retailer for information concerning the correct disposal of their old appliance.

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